

Compact 2-channel instrument device dedicated for small laboratories to perform coagulation screening and DDimers tests

Smart way to test

- 2 independent measuring channels
- Easy to operate, innovative user interface
- QC management
- Master calibration curve (no calibration needed)

Comprehensive security

- Bar code reading capability
- External data storage (USB device)
- Reagent ready-to-use

Increase efficiency

- 3 reference measuring methods : clotting, turbidimetric and chromogenic
- 120 sec maximum measurement
- Maintenance free

User friendly features

- 1 000 data storage capacity
- External printer connectivity
- Small footprint
- Multiple connectivities
 - Laboratory information system
 - Optional Yumizen G data software



Yumizen

G200

By **HORIBA**
Medical

Coagulation Yumizen Range



Yumizen G200

Dimensions & Weight:

Length : 102 mm, Height : 50, Depth : 165 mm, Weight : 0,4 kg

Throughput:

120 sec maximum measurement

Temperature:

15°C-30°C

Power:

100-240VAC, 47-63Hz, 25W
External power supply

Export:

2 x USB plug :
1 x Master USB connector : to connect an external printer, barcode reader or external storage device to the instrument.
1 x Slave USB connector : to connect the instrument to a computer as an external device.

Screen:

Color LCD touch screen 53 x 71 mm

Storage capacity:

1000 results

Parameter :

Clotting tests (640 nm)
Screening tests (PT, APTT, Fibrinogen, TT)
Factors (II, V, VII, X, VIII, IX, XI, XII)
Turbidimetric test (570 nm)
D-Dimer
Chromogenic test (405 nm)
ATIII

Consumables (kit) :

Yumizen G200	1 pc.
Power supply unit	1 pc
Magnetic mixer (2 pcs / set)	1 set
User Manual	1 pc

Norms:

EN 61326-2-6 : 2006
IEC 61010-1 : 2001
IEC 61010-2-101 : 2002
EN ISO 13485 : 2003 and 2012
NF EN ISO 14971 - 2013

Directives:

98/79/EC (DIV)
2011/65/EU (RoHS)
1907/2006 (REACH)
2014/30/EU (Electro Magnetic Compatibility)
2014/35/EU (Voltage limits)

HORIBA Medical


FRANCE +33 (0)4 67 14 15 15 - ITALY +39 / 06 51 59 22 1 - SPAIN +34 / 91- 353 30 10 - PORTUGAL +351 / 2 14 72 17 70 - UK +44 (0) 1604 542650
POLAND +48 / 22 6732022 - USA +1 / 949 453 0500 - BRAZIL +55 / 11 2923-5439 - THAILAND +66 / 2 861 59 95 - CHINA +86 / 21 3222 1818
INDIA +91 / 11 4646 5000 - GERMANY AXON LAB AG +49 / 7153 92260 - DISTRIBUTORS NETWORK +33 (0)4 67 14 15 16

HORIBA Medical online : <http://www.horiba.com/medical>

Certified
ISO 9001 / 14001
OHSA 18001
ISO 13485

FR 09/2080
FR 09/0543
FR 09/0219
FR 05/8448

Operating IMS

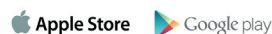

Explore the future

Automotive Test Systems | Process & Environmental | Medical | Semiconductor | Scientific

HORIBA