

BS-480

Chemistry Analyzer

- Discrete, random access, fully automated
- Constant 400 tests per hour, up to 560 tests per hour with ISE
- Up to 78 on-board chemistries & 3 ions
- Refrigerated reagent compartment at 2~10 °C
- 90 sample positions and 90 cuvettes
- Automatic probe cleaning, liquid level detection & collision protection
- 8-step auto washing system
- 12 wavelengths: 340~800nm
- Minimum reaction volume 120 μ l
- Intelligent auto-recovery from probe collision
- Advanced multi-functional and easy-to-use software
- Enhanced liquid system with easy maintenance
- Dedicated system pack reagents



Experience the Added Value

mindray
healthcare within reach

BS-480

Chemistry Analyzer

Mindray solution for clinical chemistry laboratories

Enhanced Performance

Minimum sample volume: 1.5 μ l

Reagent volume: 10~350 μ l


Minimum reaction volume: 120 μ l

Probe carry-over: $\leq 0.05\%$

Intelligent auto-recovery from probe collision without system stopping

Absorbance range: 0~3.3 A (10mm conversion)

System auto-start function

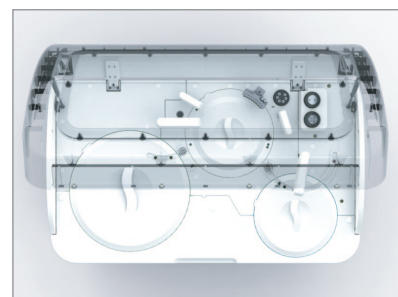

Advanced Software

Inherits the latest software platform and UI design

Supports touch screen interface

Multi-functional and easy to use

Auto re-run function

Auto serum index

Real-time QC status monitoring

Test summary function

Reflex function



Dedicated System Pack Reagents

48 parameters of dedicated reagents available

More than 17 others are coming

Original calibrators with metrological traceability for all parameters

Original controls for all parameters



Easy Maintenance

Optimized liquid system

All containers & kits are located in the front of the machine

Easy for replacing, maintenance and trouble shooting

Step-by-step maintenance guide in the software



DISTRIBUTOR:



MINDRAY is a trademark of Shenzhen Mindray Bio-Medical Electronics Co., Ltd. Specifications subject to changes without prior notice.
© 2012 Shenzhen Mindray Bio-Medical Electronics Co., Ltd. All rights reserved.
P/N: ENG-BS480-210285x2-20120327

mindray

Mindray is listed on the NYSE under the symbol "MR"
Mindray Building, Keji 12th Road South, High-tech Industrial Park,
Nanshan, Shenzhen 518057, P.R. China
Tel: +86 755 81888998 Fax: +86 755 26582680
E-mail: intl-market@mindray.com Website: www.mindray.com